

201603229001

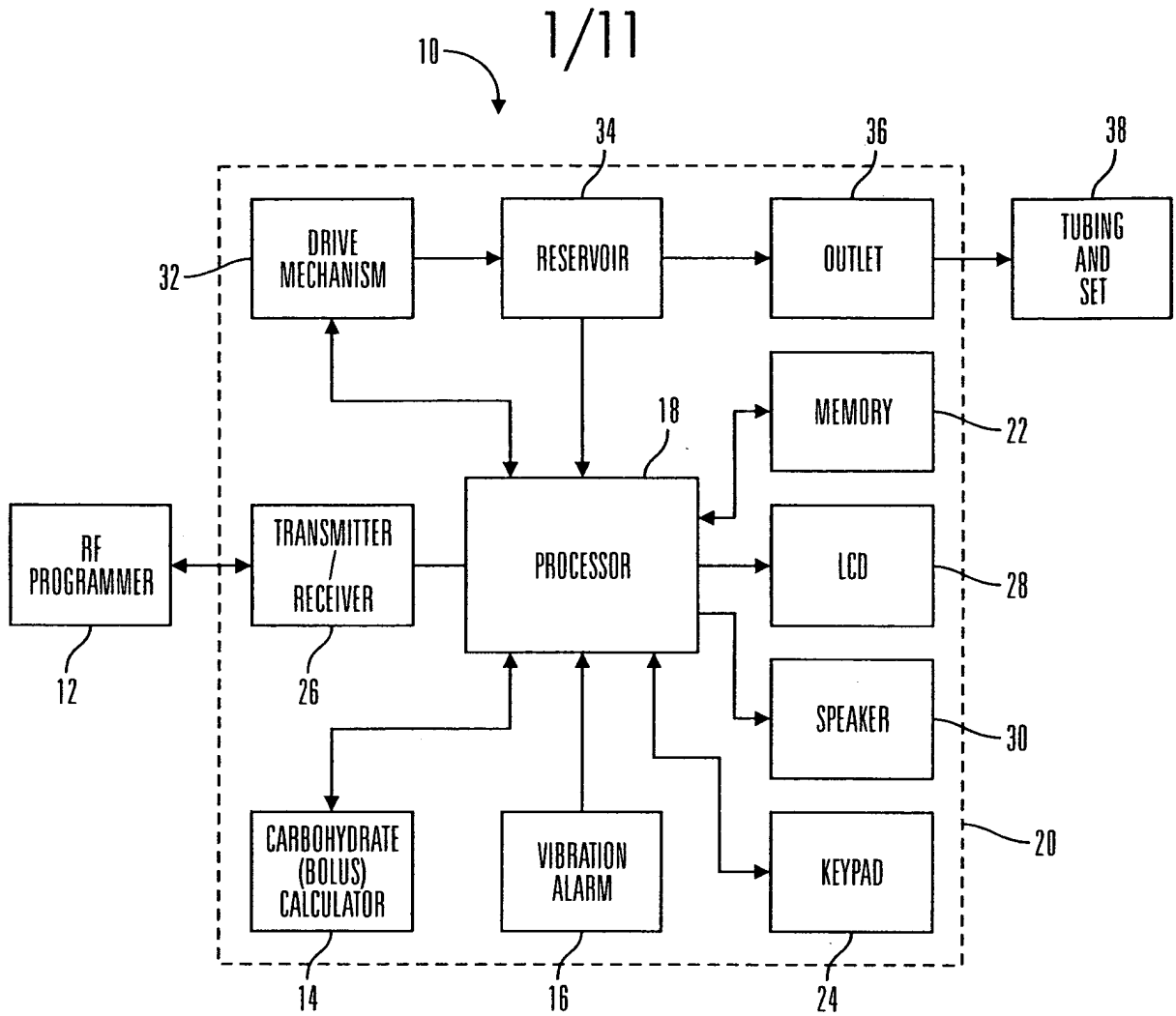


FIG. 1

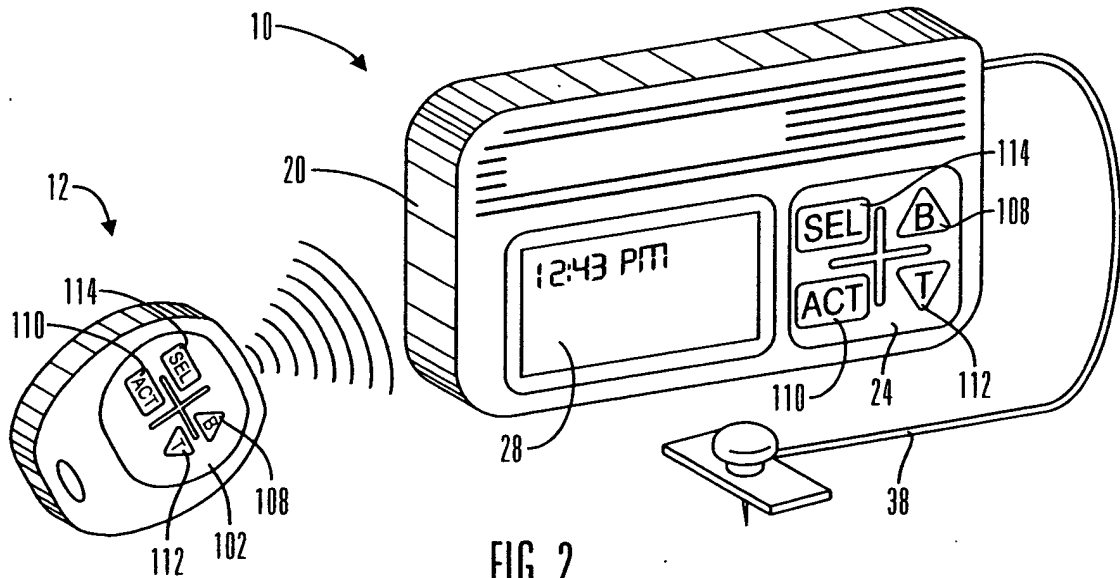


FIG. 2

2/11

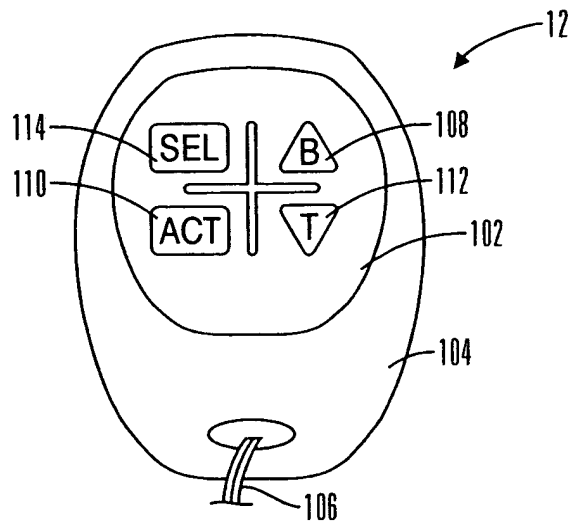


FIG. 3

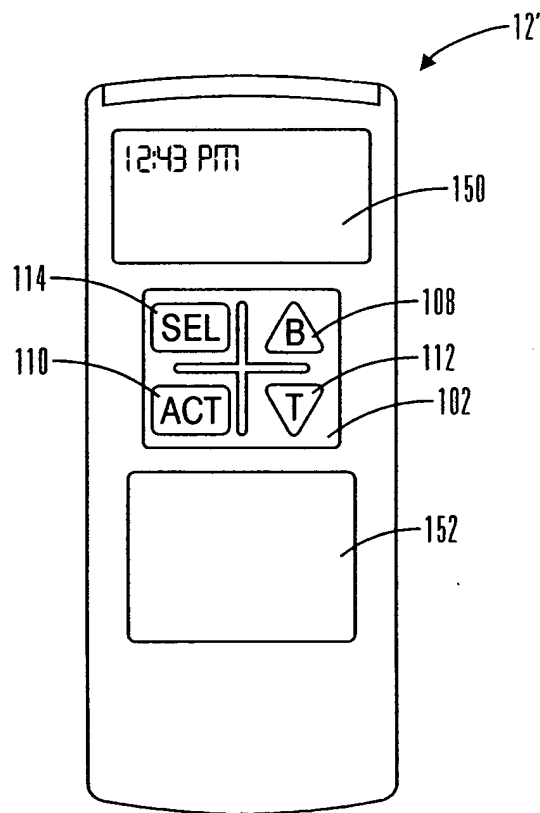


FIG. 4

3/11

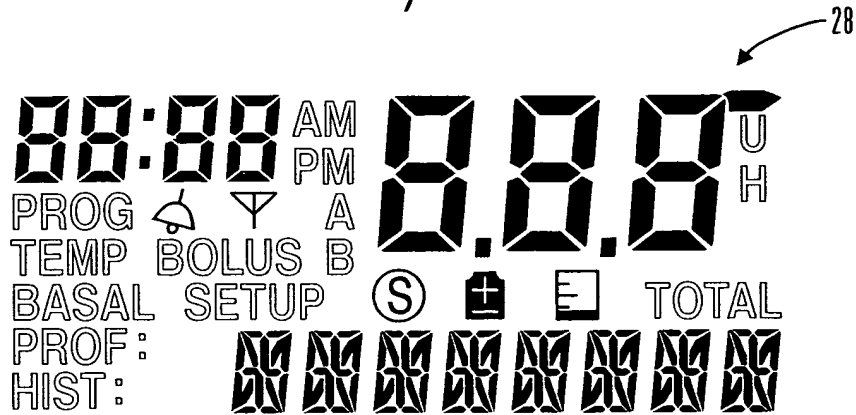


FIG. 5

NO.	SETUP II MENU OPTIONS
1	VARIABLE BOLUS MODE ON/OFF SCREEN
2	BOLUS ESTIMATOR ON/OFF SCREEN
3	MAXIMUM BOLUS SCREEN
4	MAXIMUM BASAL SCREEN
5	ALARM REVIEW SCREEN
6	RF DEVICE ID SCREEN
7	AUDIO BOLUS MODE ON/OFF SCREEN
8	PERSONAL DELIVERY PATTERN SCREEN
9	ALARM MODE (AUDIO/VIBRATE) SCREEN
10	CHILD BLOCK ON/OFF SCREEN
11	TIME DISPLAY SCREEN
12	INSULIN CONCENTRATION SCREEN
13	SETUP I SCREEN
14	SETUP EXIT SCREEN

FIG. 6

10062838 013102

508 PUMP - BOLUS SEQUENCE MENU

4/11

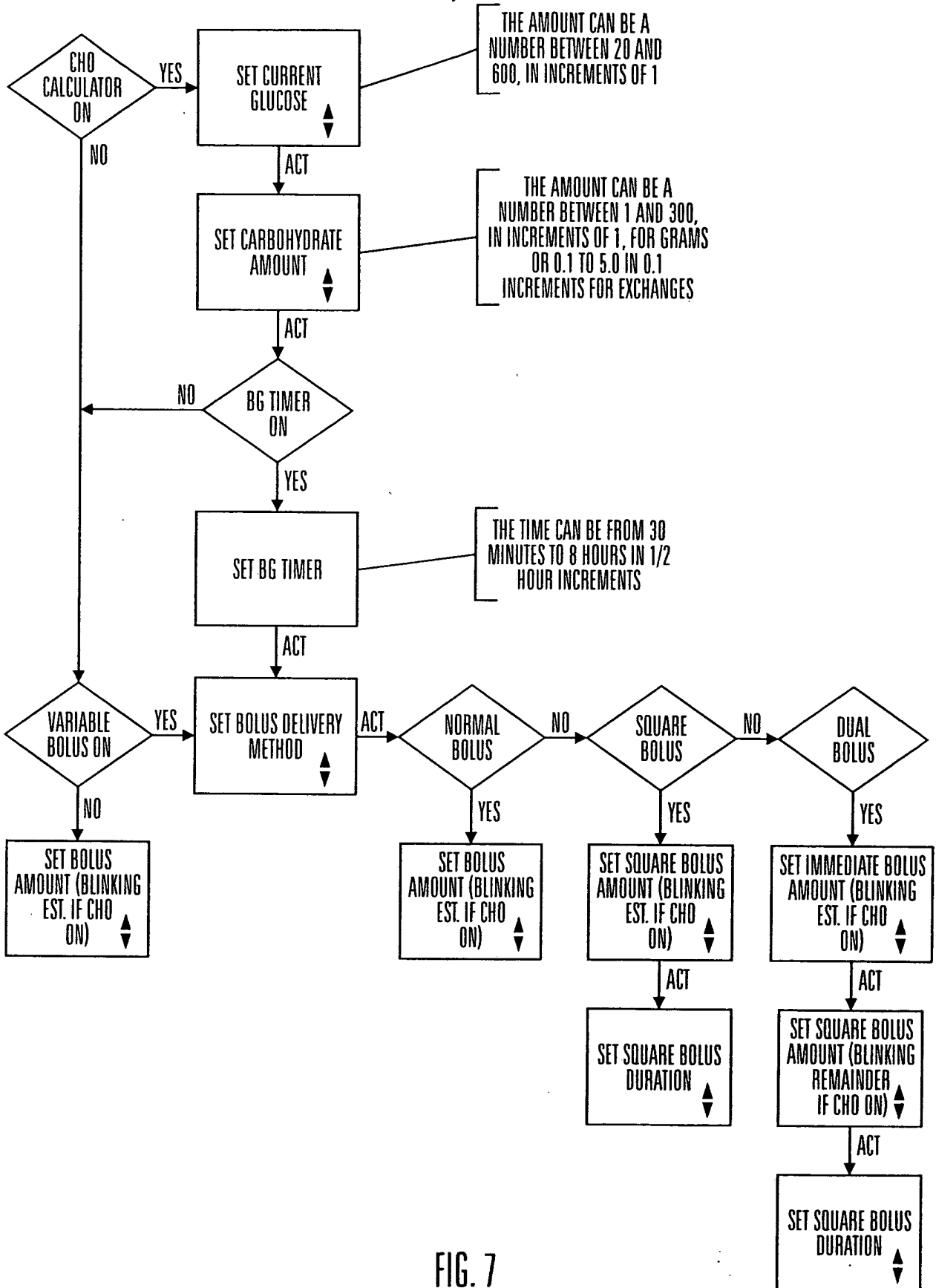


FIG. 7

5/11

508 PUMP - SETUP II MENU

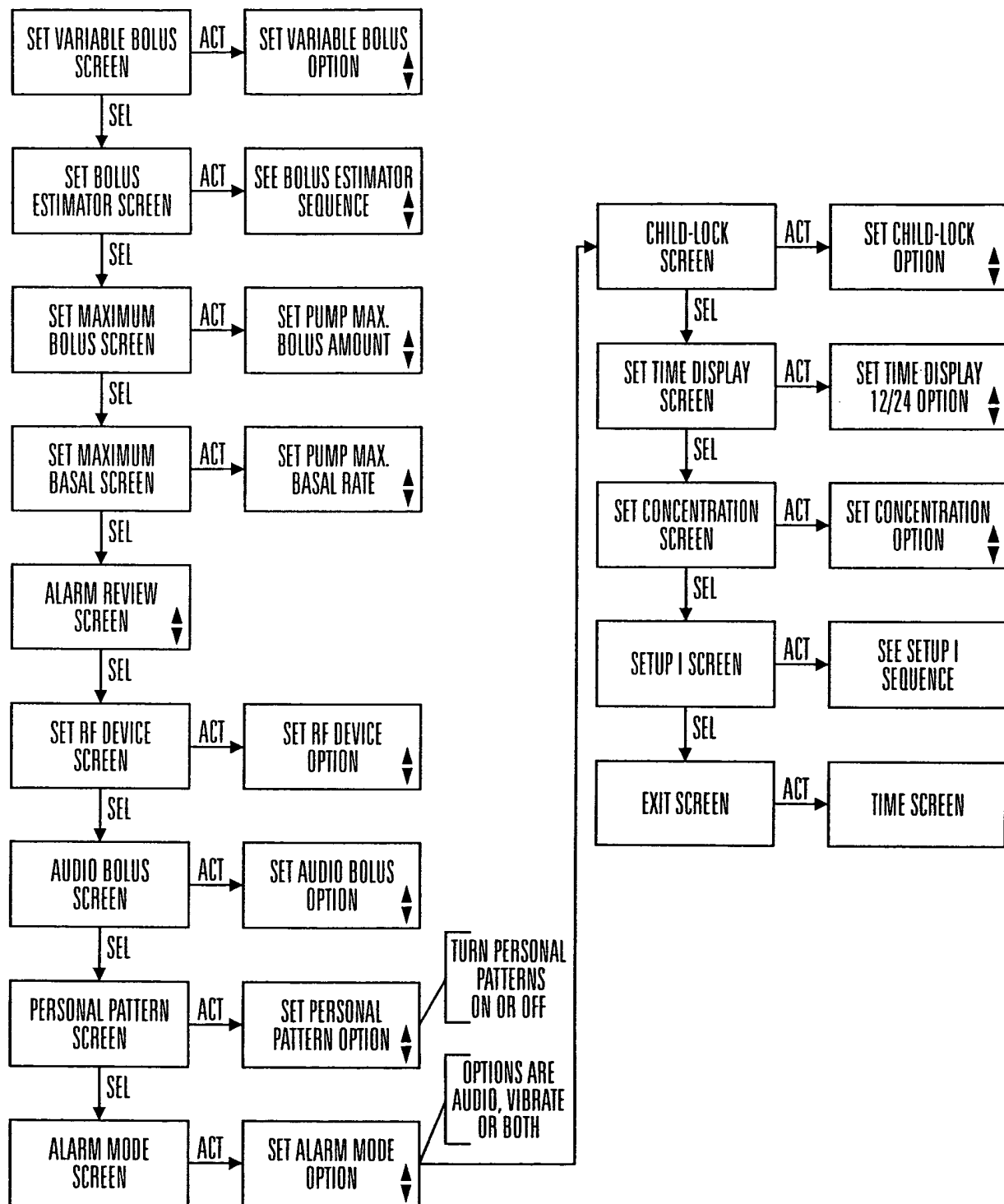


FIG. 8(a)

6/11

508 PUMP - BOLUS ESTIMATOR SEQUENCE

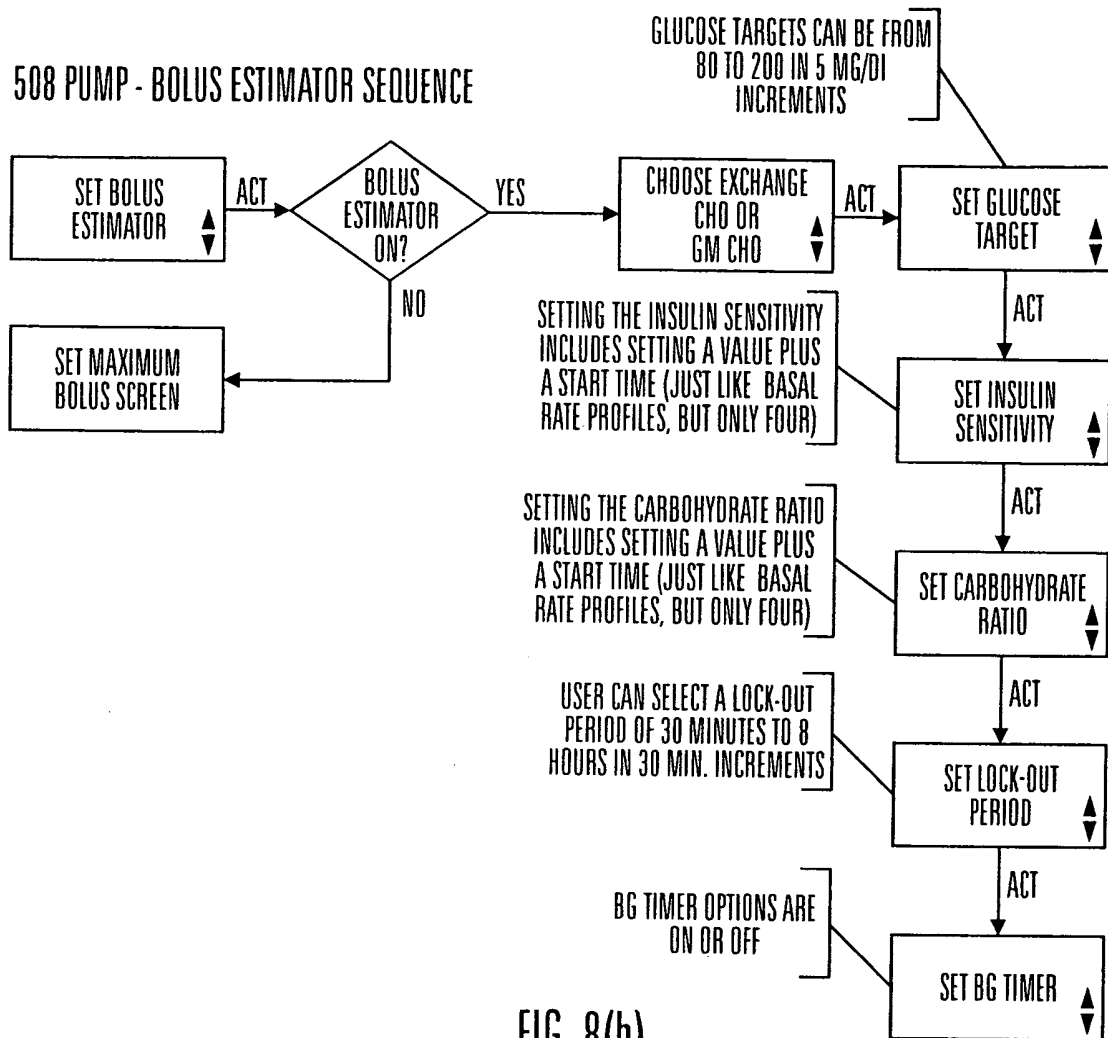


FIG. 8(b)

NO.	MAIN MENU OPTIONS
1	TIME SCREEN
2	BOLUS & BOLUS HISTORY SCREEN
3	SUSPEND/RESUME PUMP SCREEN
4	BASAL RATE SCREEN
5	TEMP. BASAL RATE SCREEN
6	DAILY TOTALS REVIEW SCREEN
7	PRIME BOLUS SCREEN
8	SETUP I SCREEN
9	SETUP II SCREEN

FIG. 9

7/11

NO.	SETUP I MENU OPTIONS
1	TIME & DATE ADJUSTMENT SCREEN
2	AUTO-OFF SCREEN
3	BEEP VOLUME SCREEN
4	SELF-TEST SCREEN
5	SETUP II SCREEN
6	SETUP EXIT SCREEN

FIG. 10

508 PUMP - MAIN MENU

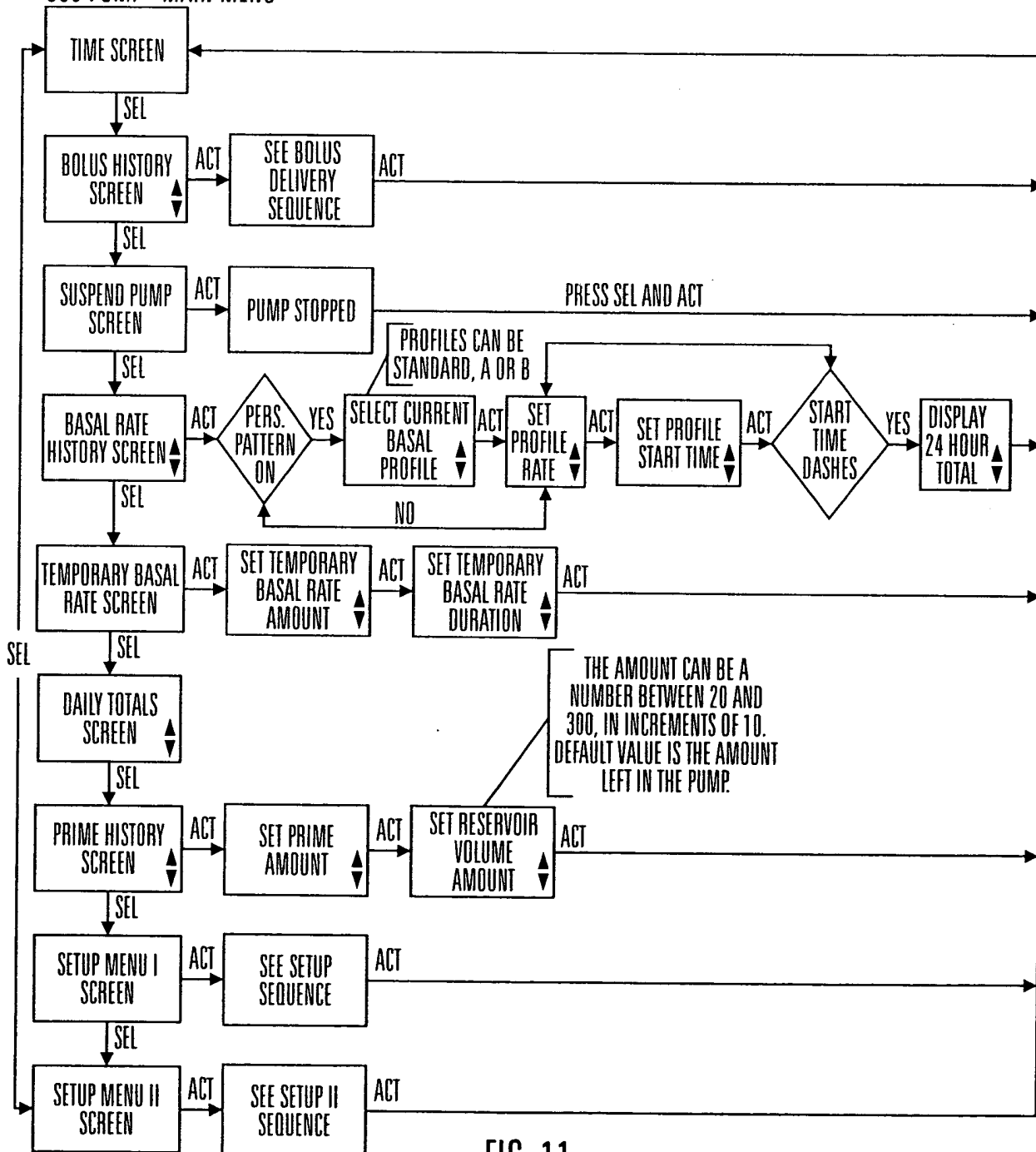


FIG. 11

+

8/11

508 PUMP - SETUP I MENU

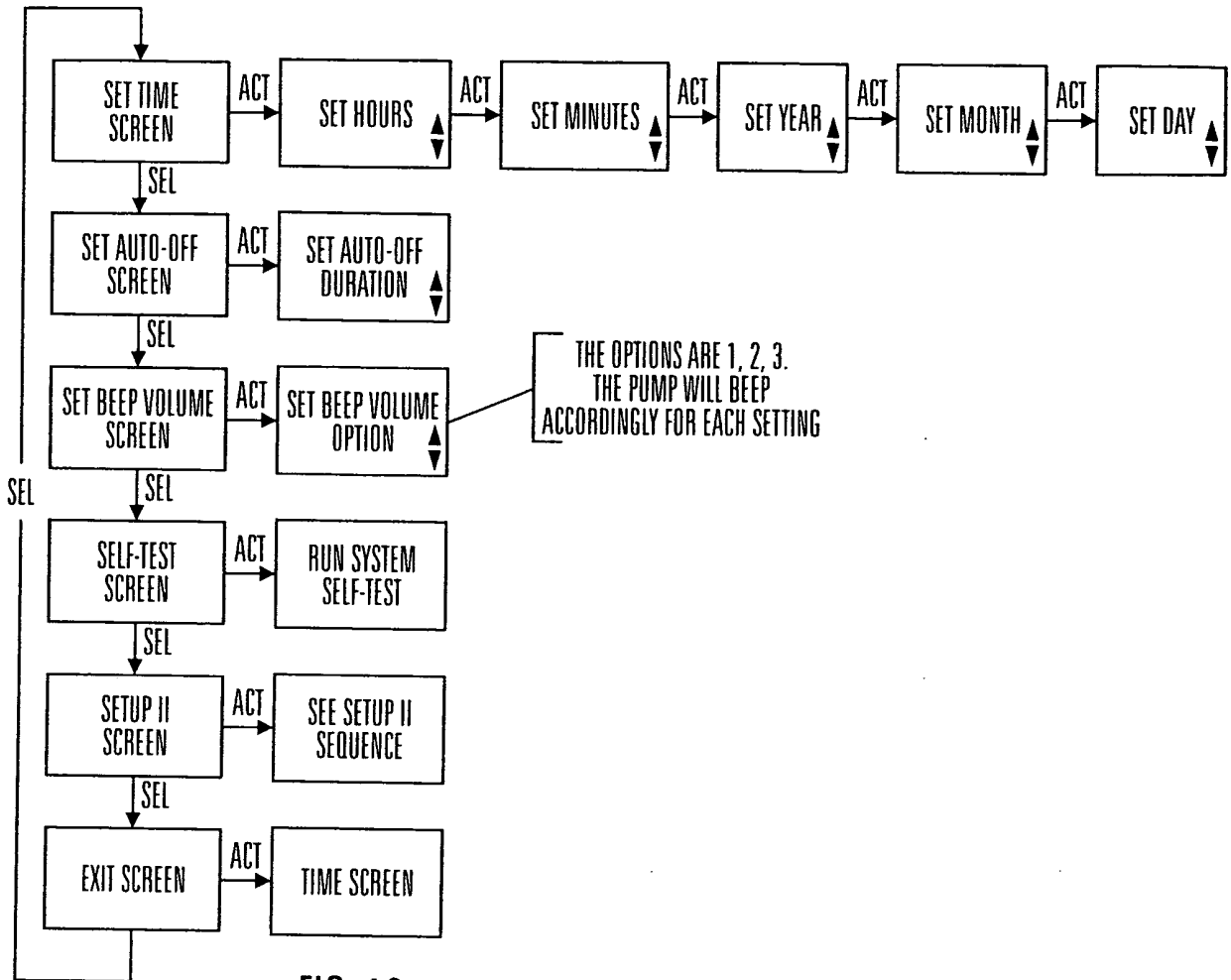


FIG. 12

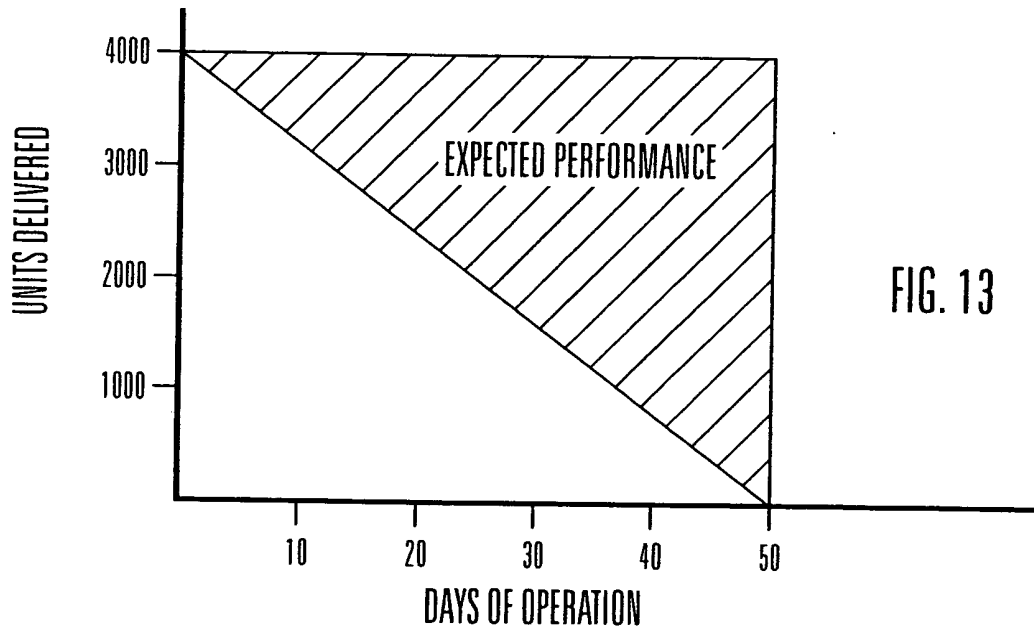


FIG. 13

2025-03-04 10:00:00

+

9/11

PARAMETER	FACTORY DEFAULT
INSULIN CONCENTRATION	U100
BASAL RATE PROFILES	PROFILE 1, 0.0 U/H (PERSONAL PATTERNS 2 AND 3 ARE ALSO SET TO PROFILE 1 OF 0.0 U/H). PERSONAL PATTERN 1 IS THE DEFAULT.
MAXIMUM MEAL BOLUS	10.0 U
MAXIMUM BASAL RATE	2.0 U/H
TIME	12:00 AM, JAN 1, 1998 AT TIME OF PUMP POWER ON RESET (E-01)
TIME DISPLAY FORMAT	12 HOUR (FOR U.S.A.), 24 HOUR (FOR EUROPE)
AUDIO BOLUS MODE	OFF
AUDIO BOLUS INCREMENT	0.5 U (U100)
BEEP VOLUME	2 (MEDIUM)
AUTOMATIC OFF DURATION	-- (DISABLED)
DAILY TOTALS	ALL 90 DAYS 0.0 U (USER CAN REVIEW LAST 7 ONLY)
BOLUS HISTORY	ALL 450 BOLUSES 0.0 U (USER CAN REVIEW LAST 24 ONLY)
PRIME BOLUS HISTORY	ALL 50 BOLUSES 0.0 U (USER CAN REVIEW LAST 9)
ALARM MODE	AUDIO
LOW VOLUME ALERT	OFF
ALARM CLOCK	OFF
CHILD LOCK	OFF
BOLUS ESTIMATOR	OFF
TARGET BLOOD GLUCOSE	100 MG/DL
INSULIN SENSITIVITY	30
CARBOHYDRATE RATIO	15
EXCHANGE RATIO	1.0
BOLUS ESTIMATOR LOCKOUT	1 HR
BLOOD GLUCOSE TIMER	OFF
RF-DEVICE	OFF
ALARM HISTORY	BLANKS (USER CAN REVIEW LAST 12 ONLY)

FIG. 14

10062932 "043402

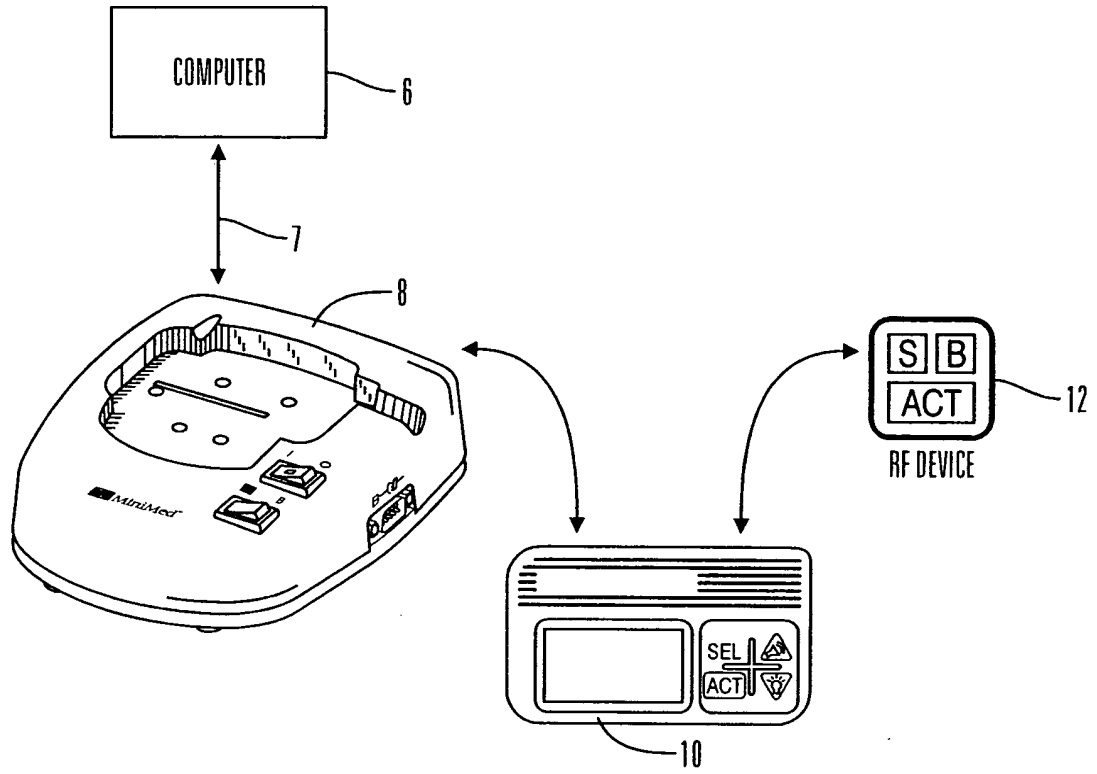


FIG. 15

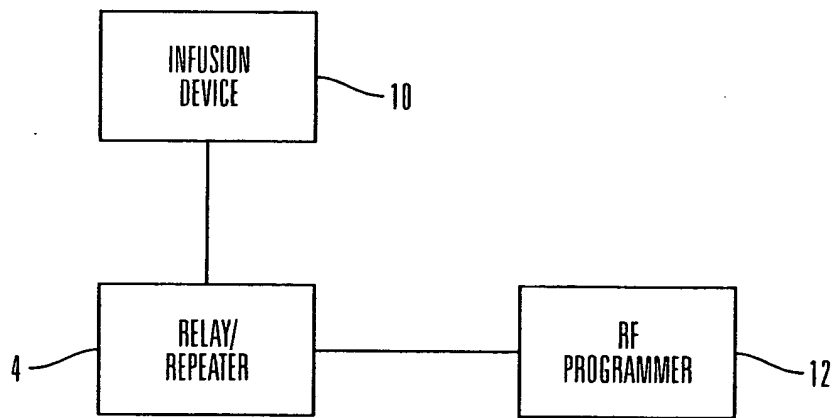


FIG. 16

10062838 013102

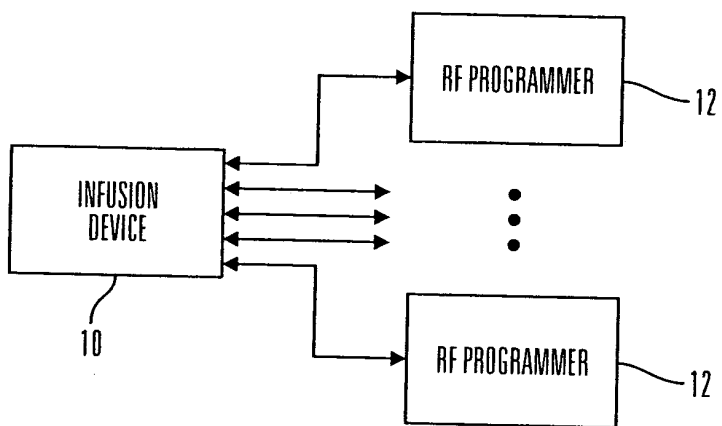


FIG. 17